

EDGEONYX™

Take The EdgeOnyx Challenge

Push it. Stress it. Test its limits. We dare you.



EDGEENDO®

Official AAE Launch

EDGEONYX™



The EdgeOnyx rotary file system features a minimally invasive design (<0.8 max flute diameter) and a simplified sequence that guides clinicians through each stage of root canal preparation.

EdgeOnyx is based on the time-tested design of the original EdgeFile®, but with several design modifications to improve performance such as, advanced proprietary heat treatment to maintain flute cutting edge fidelity.

The most notable upgrade is a stiffened shank portion, which reduces “wobble” and improves visibility and tactile control during treatment. Powered by a proprietary heat treatment process, EdgeOnyx delivers maximum flexibility without sacrificing shank stiffness. The result is enhanced control, smoother instrumentation and performance you can rely on case after case.

Another significant improvement is the anti-ledging safety tip which helps the files navigate even the most challenging canals. The black shank makes the calibration marks more noticeable for improved depth control.


The system protocol has been streamlined to guide the clinician through a “find, follow, finish” approach.

The system features a unique 15mm #09.09 orifice shaper (X-Find), size #15 or 17.04 glidepath file options (X-Follow) and is rounded out with #20-45.04 (X-Finish) files.



A parabolic cross-section delivers efficient cutting and enhanced debris removal.

Simplified Technique

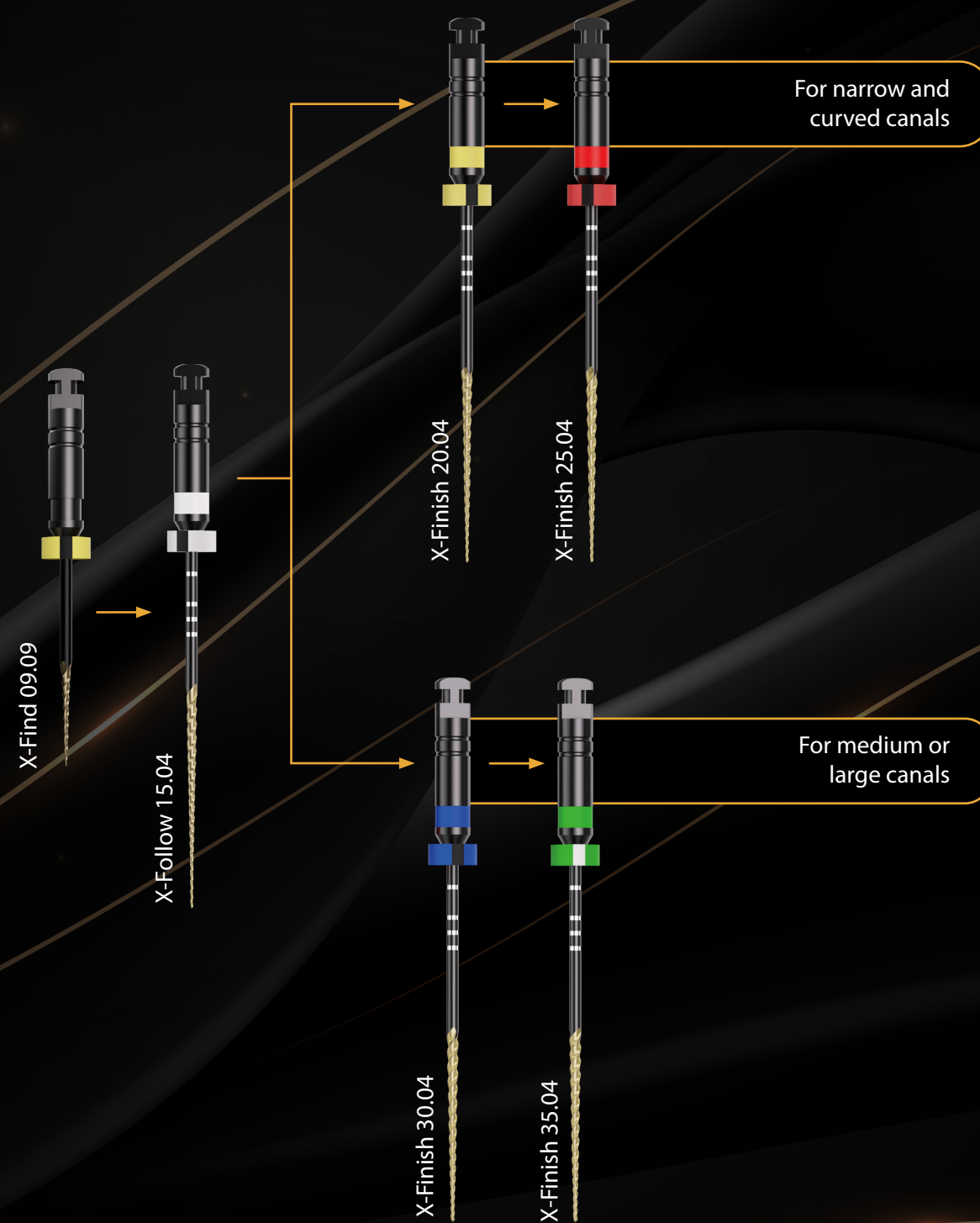

Triflex System



X-Find (Orifice Opener)

Designed for precise, minimally invasive endodontics. Conical design of the X-Find maximizes coronal dentin preservation



X-Follow (Glidepath)

Optimizes the glidepath trajectory for a smoother transition, enhanced efficiency and safety, and minimized risk of ledge formation.

Available in two size options: #15.04 or #17.04



X-Finish (Shaping)


Range of sizes allowing the completion of preparation in different clinical cases.

Available in tip sizes ranging #20-45.04 curated to meet the needs of most clinical scenarios.

EDGEONYX™

\$28.80* / 6pk

Balanced Performance

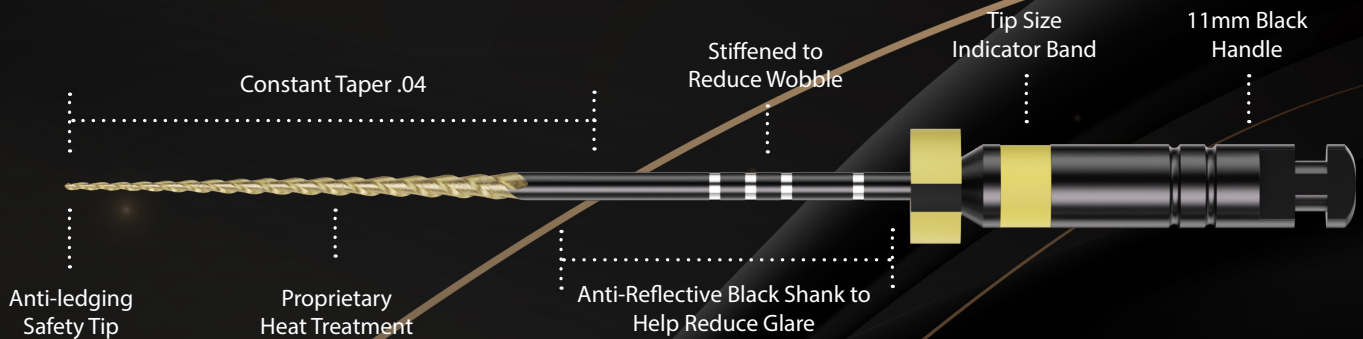

- Parabolic Cross Section
- Available in 2 Lengths: 21 & 25mm w/ 15mm orifice opener**
- Speed: 300-500 RPM
- Torque: 2.5 Ncm



Precise Control



Pre-Bendable



Minimally Invasive

X-FIND Orifice Opener - 6 Pack

Size	Color	15mm only
09.09	YELLOW	EO090915HTX

6 Pack			
Size	Color	21mm	25mm
X-FOLLOW Glidepath Options			
15.04	WHITE	EO150421HTX	EO150425HTX
17.04	PURPLE	EO170421HTX	EO170425HTX

X-FINISH Shaping Options

20.04	YELLOW	EO200421HTX	EO200425HTX
25.04	RED	EO250421HTX	EO250425HTX
30.04	BLUE	EO300421HTX	EO300425HTX
35.04	GREEN	EO350421HTX	EO350425HTX
40.04	BLACK	EO400421HTX	EO400425HTX
45.04	WHITE	EO450421HTX	EO450425HTX

Assorted Packs - \$28.80/ 6 Pack

Size	Assorted Packs - \$28.80/ 6 Pack	
	21mm	25mm
X-Find (20), X-Follow (15), X-Finish (20, 25, 30)	EOASST21HTX	EOASST25HTX

*Not Available For Discounts **15mm orifice opener available in #09.09



EDGEENDO®

Scan to register for Edge's Apex Academy to try EdgeOnyx

